


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Hydroxychloroquine Sulfate

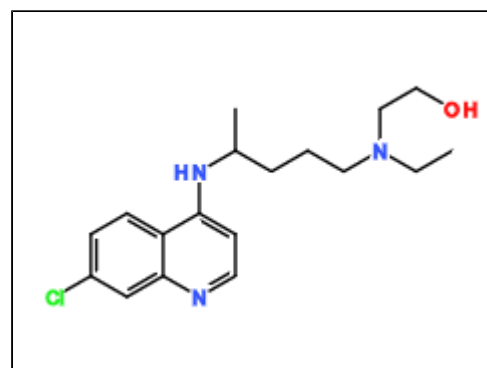

## CHROMATOGRAM

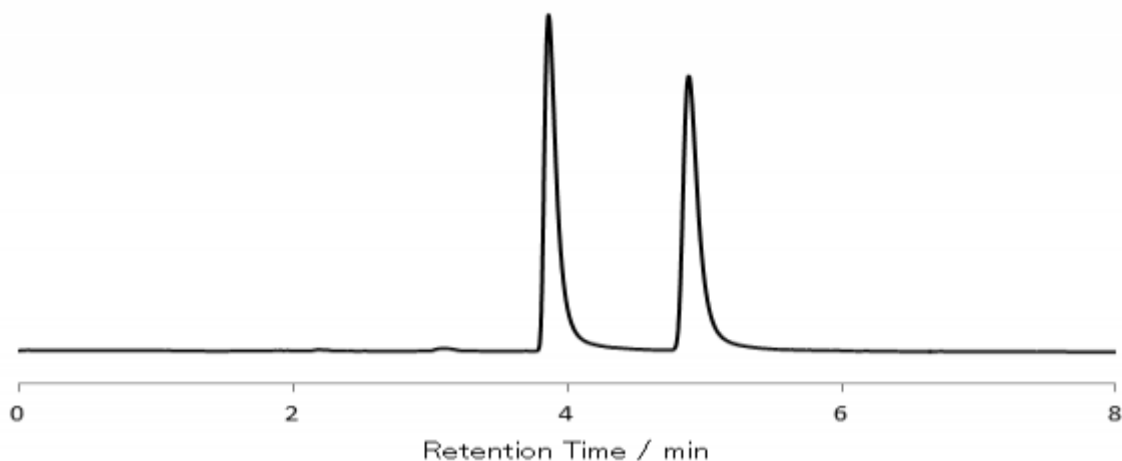

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IG-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	acetonitrile / methanol / ammonium acetate (weight) = 900 / 100 / 0.33		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)	0.003 (mg)	Sample Concentration	
Detection	UV-VIS 328 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.9	1.10	1 - 2	5.68
2	4.9	1.70		