


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

2,2'-Bis(bromomethyl)-1,1'-binaphthalene

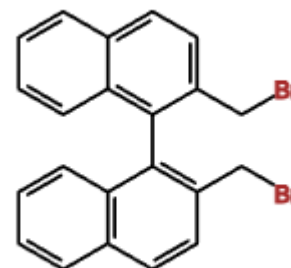

## CHROMATOGRAM

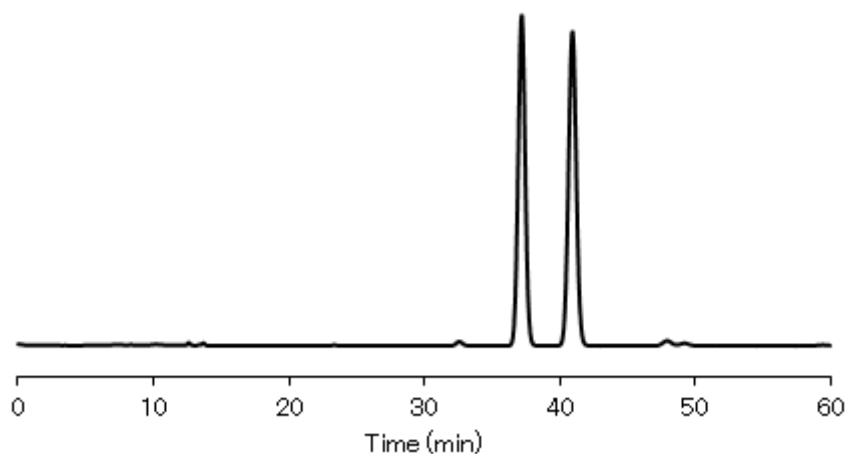

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IH-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	water / acetonitrile = 40 / 60		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	10	Sample Concentration	0.5 mg/mL
Detection	UV-VIS 285 nm		
Elution Order	第一ピーク : R体 第二ピーク : S体		

No.	Rt	K'	Alpha Value	Resolution
1	37.2		1 - 2	3.69
2	41.0			