


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ibuprofen

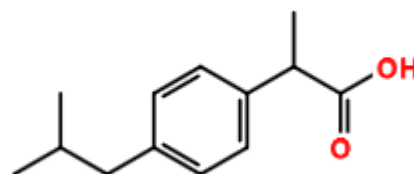

CHROMATOGRAM

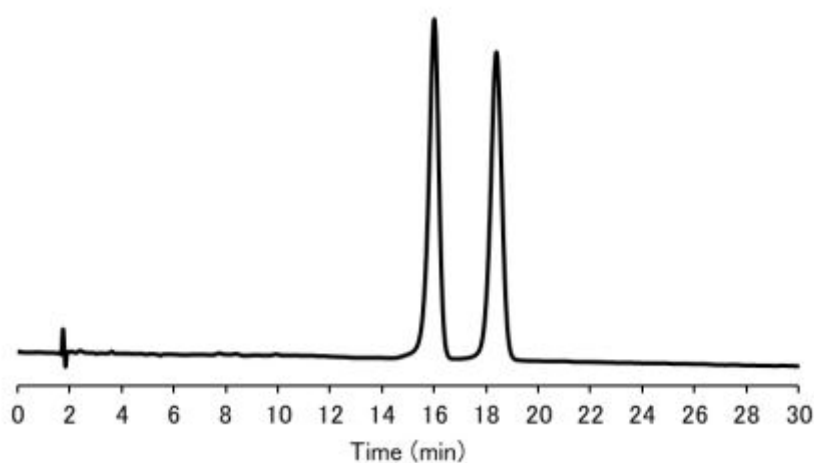

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IG-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	formic acid aq. pH2.0 / acetonitrile = 65 / 35		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)	10	Sample Concentration	1 mg/mL
Detection	UV-VIS 265 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	16.0		1 - 2	3.21
2	18.4			