


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Thalidomide

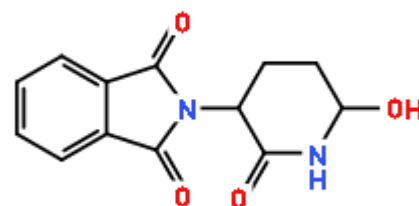

CHROMATOGRAM

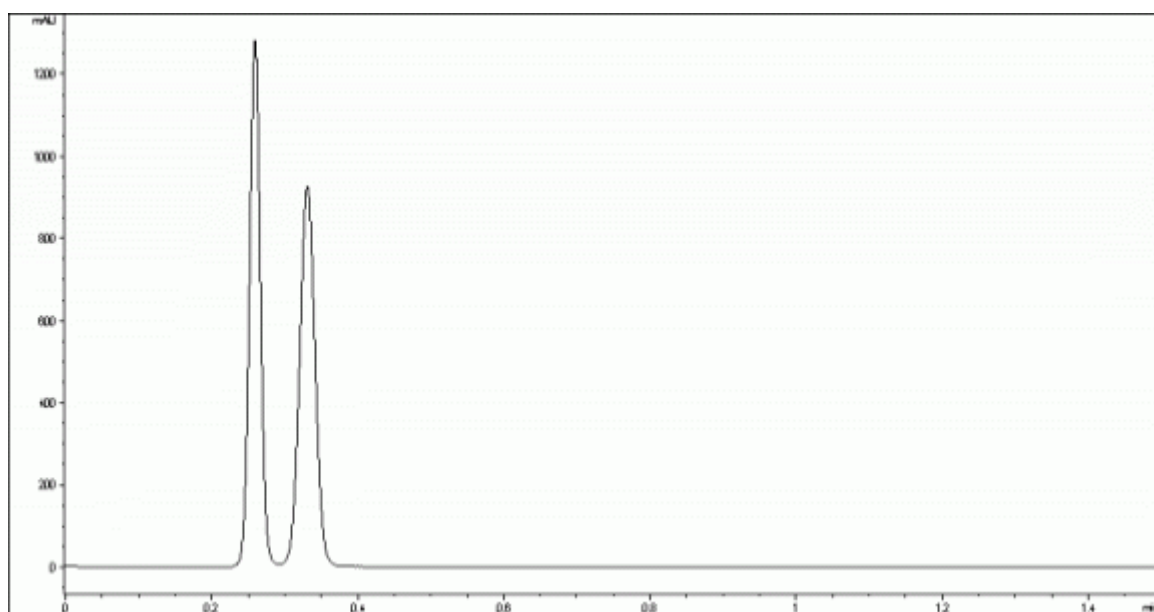

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL [®] OZ-3		
Column Size (Particle Size)	4.6 x 50 mm (3 μm)		
Mobile Phase Composition	acetonitrile = 100		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	4	Temperature (°C)	25
Injection Volume (μL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.3	0.69	1 - 2	2.19
2	0.3	1.16		