


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Clopidogrel

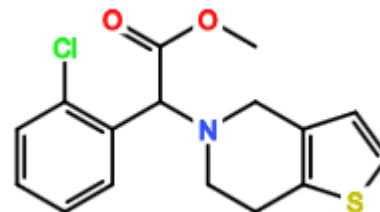

## CHROMATOGRAM

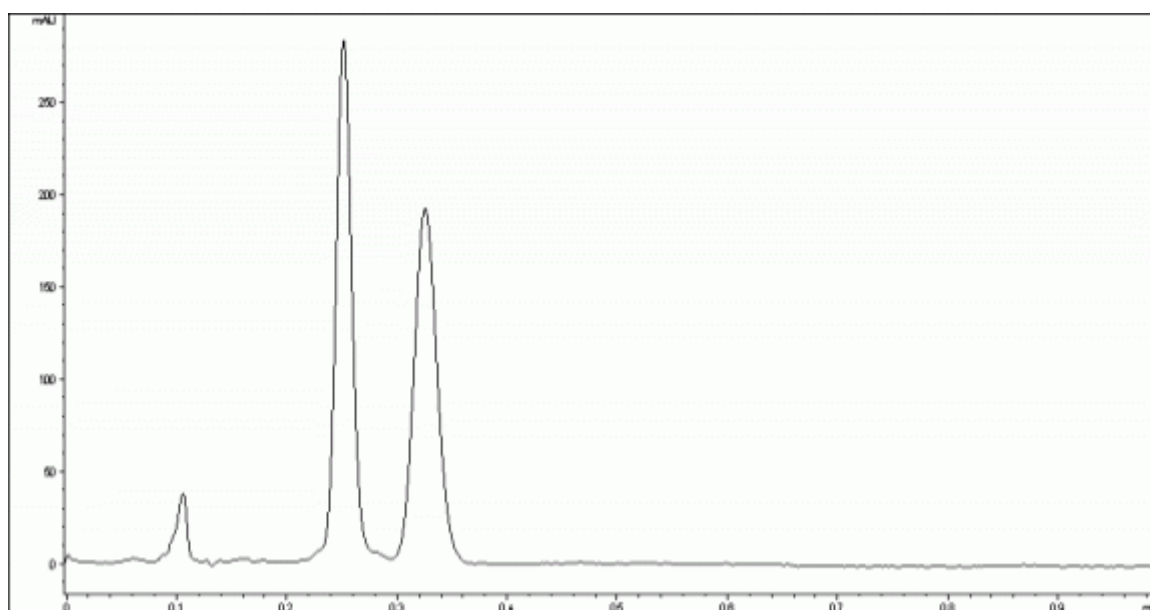

## ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OJ-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	methanol / diethylamine = 100 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	5	Temperature (°C)	25
Injection Volume (µL)	0.004 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.3	1.11	1 - 2	2.24
2	0.3	1.73		