


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

N-alpha-Benzoyl-DL-arginine-beta-naphtylamide

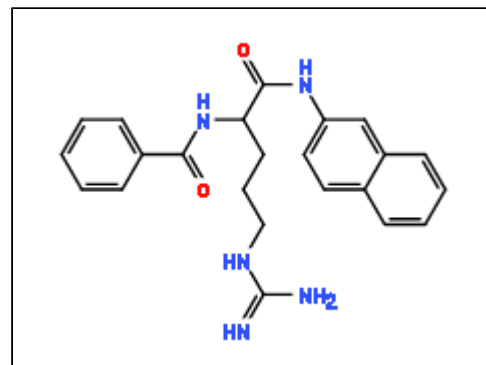

CHROMATOGRAM

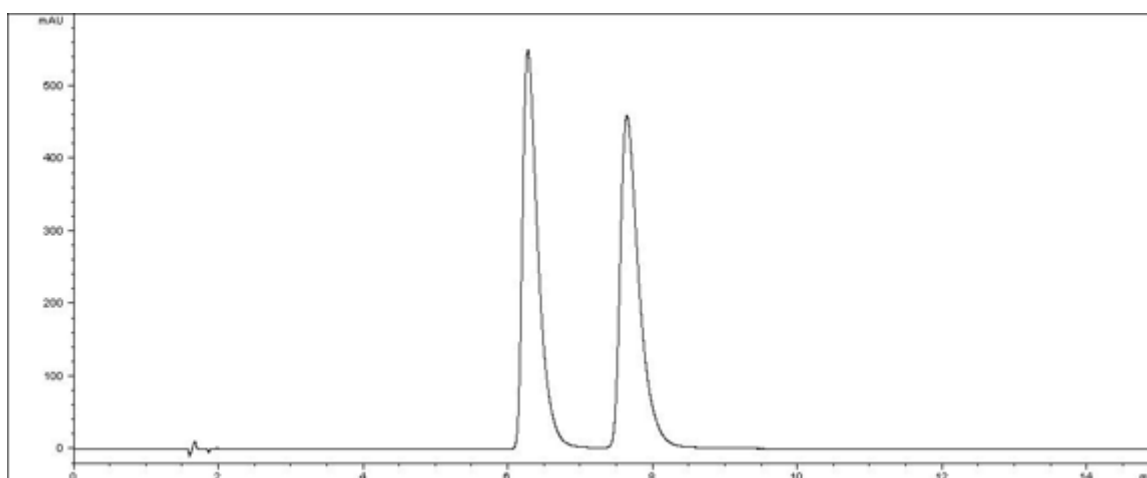

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 90 / 10		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.3	2.79	1 - 2 1.30	3.31
2	7.6	3.62		