


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Dapoxetine

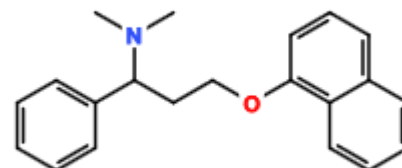

## CHROMATOGRAM

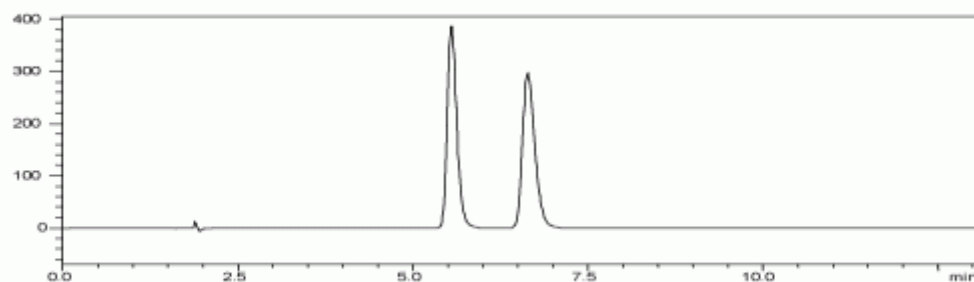

## ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OJ-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	methanol / ethanol / diethylamine = 90 / 10 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 290 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.5	2.08	1 - 2	3.52
2	6.6	2.68		