


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

gamma-(4-Fluorophenyl)-gamma-butyrolactone

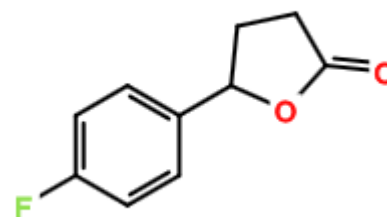

## CHROMATOGRAM

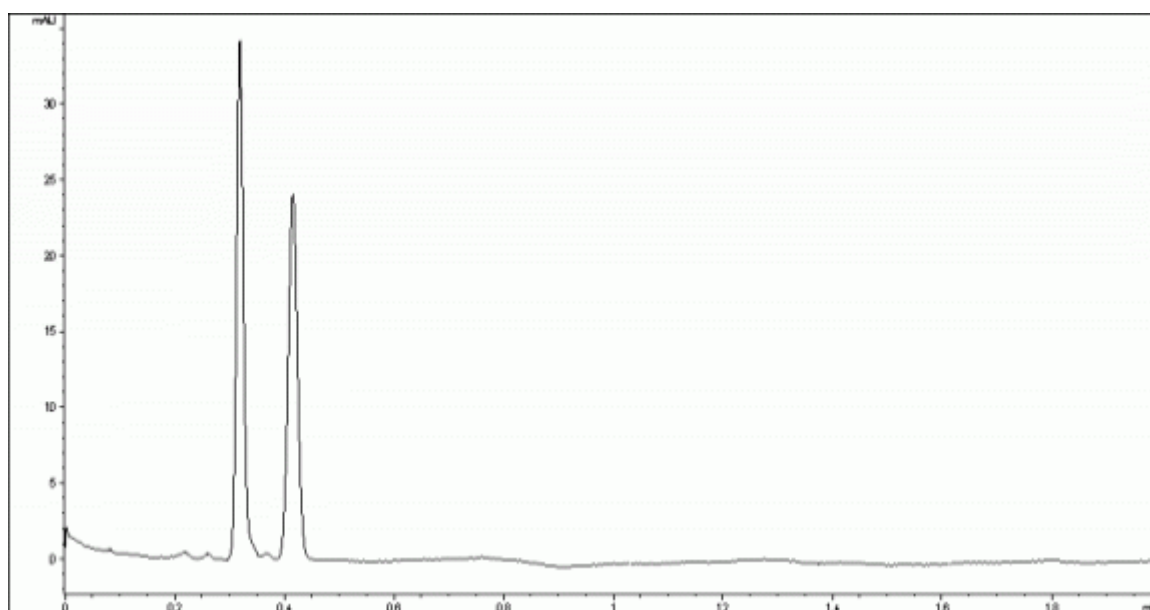

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	ethanol = 100		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	3	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.3	0.70	1 - 2	3.56
2	0.4	1.21		