


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Mephobarbital

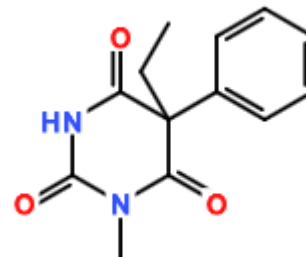

## CHROMATOGRAM

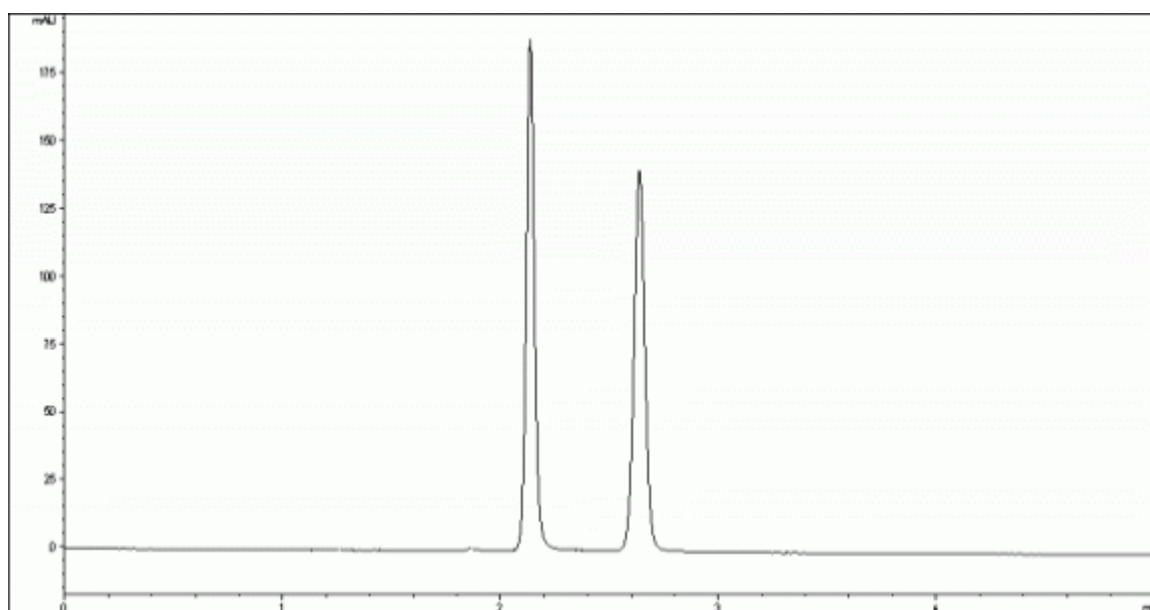

## ANALYTICAL REPORT

Column Trade Name	CHIRALCEL <sup>®</sup> OJ-3		
Column Size (Particle Size)	4.6 x 150 mm (3 $\mu$ m)		
Mobile Phase Composition	acetonitrile = 100		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature ( $^{\circ}$ C)	25
Injection Volume ( $\mu$ L)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.1	0.20	1 - 2	6.48
2	2.6	0.49		