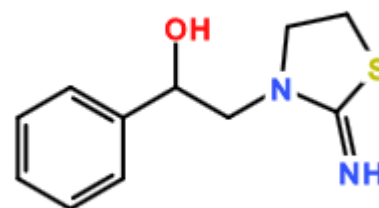



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name

2-Imino- α -phenyl-3-thiazolidine ethanol

CHROMATOGRAM

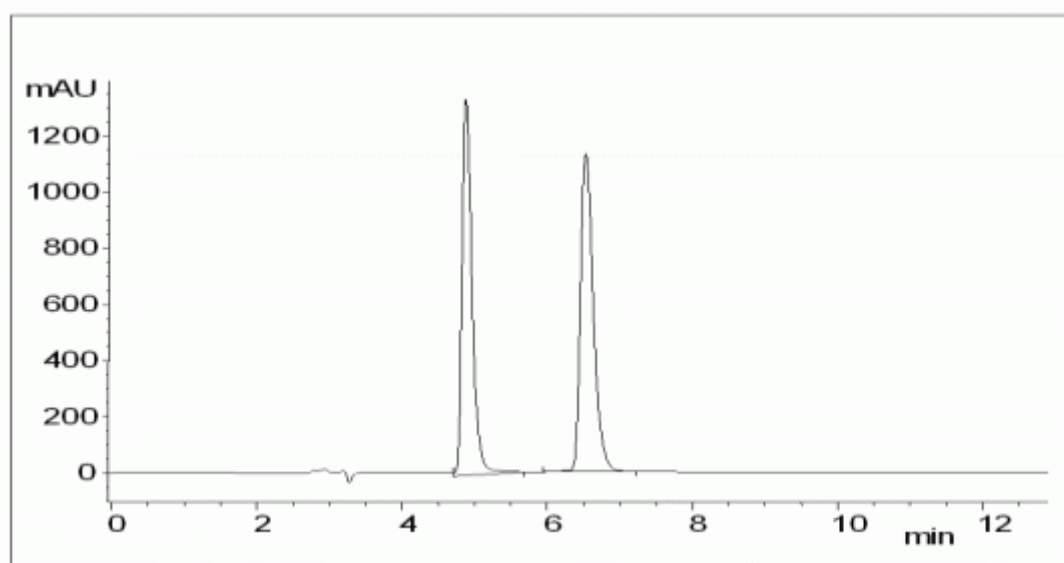

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-H		
Column Size (Particle Size)	4.6 x 250 mm (5 μ m)		
Mobile Phase Composition	ethanol / diethylamine = 100 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature ($^{\circ}$ C)	25
Injection Volume (μ L)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.9	0.63	1 - 2	5.98
2	6.5	1.18		