


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Methyl-(phenylsulfinyl)-acetate

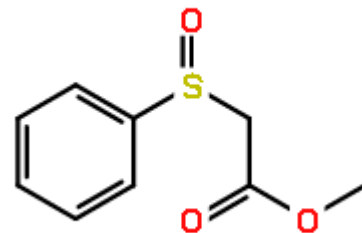

## CHROMATOGRAM

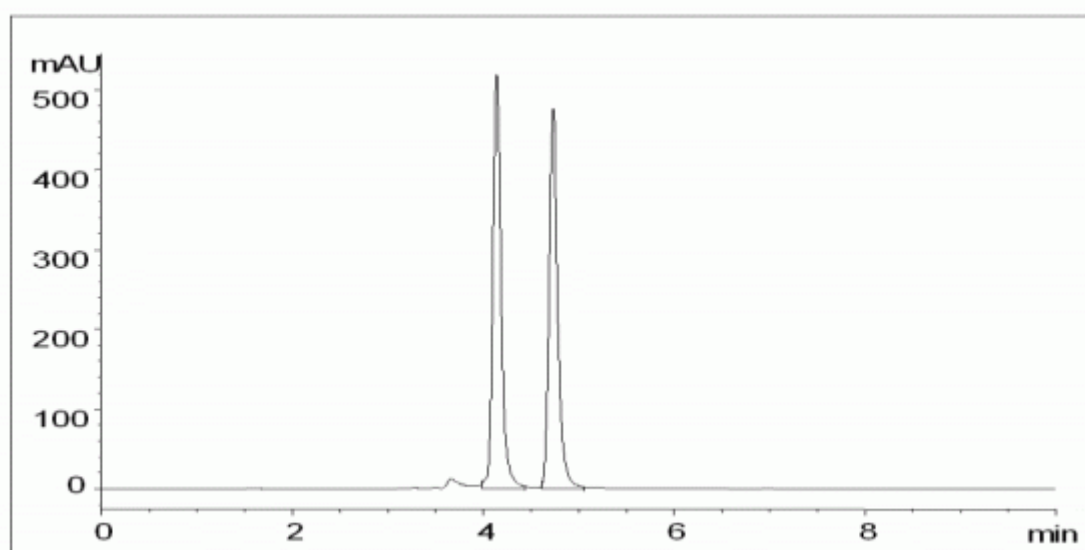

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	acetonitrile / 1-butanol = 95 / 5		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.1	0.38	1 - 2	3.60
2	4.7	0.58		