


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

N-CBZ-DL-Phenylalanine

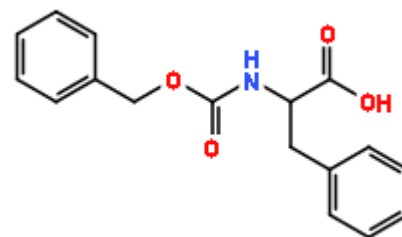

## CHROMATOGRAM

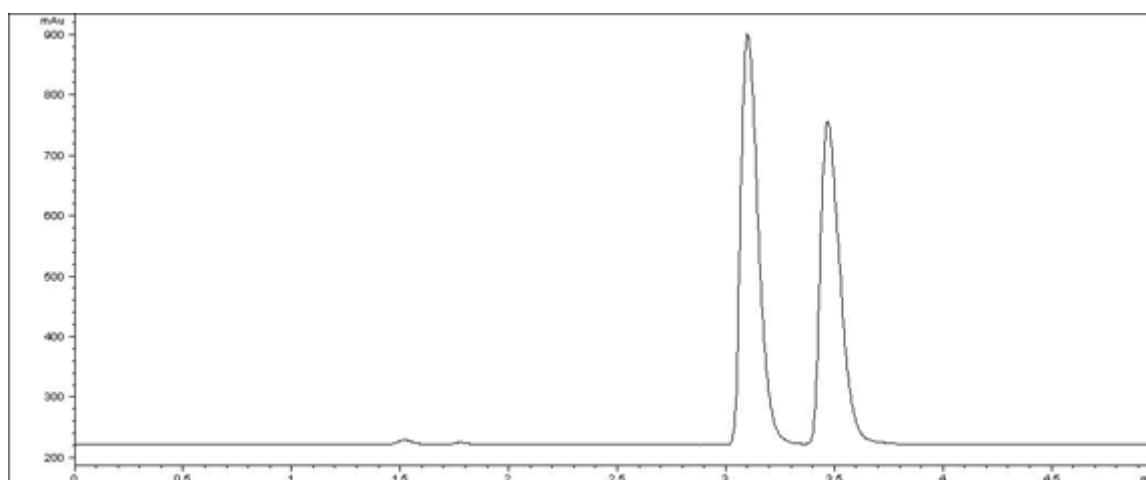

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	25mM acetic acid in methanol / water = 98 / 2		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order	D/L		

No.	Rt	K'	Alpha Value	Resolution
1	3.1	0.72	1 - 2	2.36
2	3.5	0.93		