


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

N-Acetyl-DL-Leucine

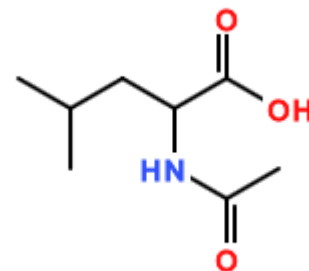

CHROMATOGRAM

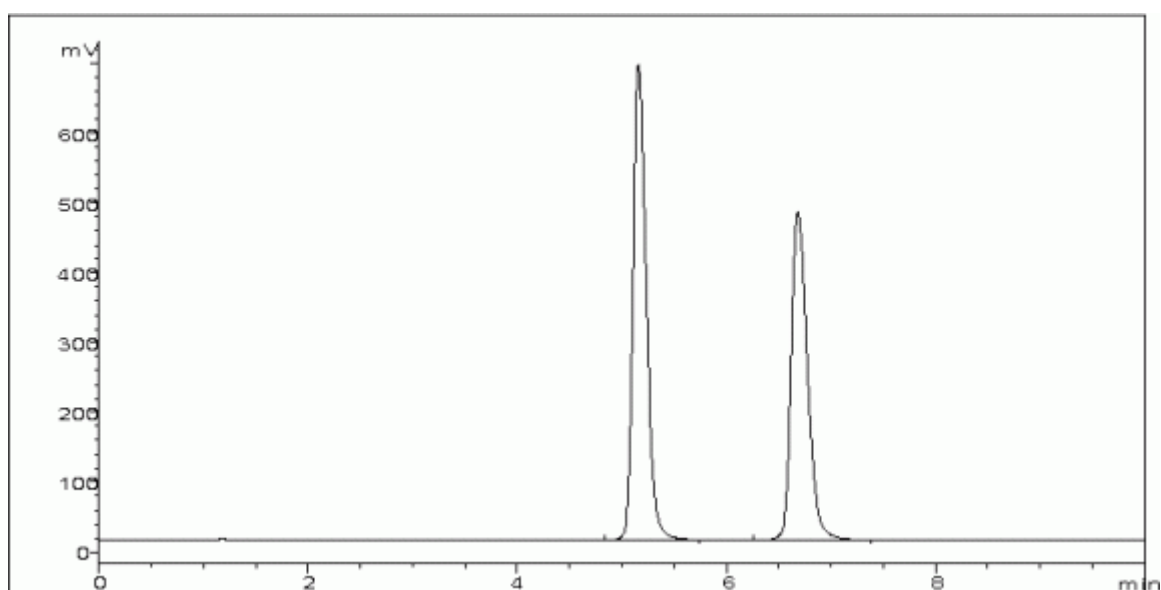

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® QD-AX		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	methanol / acetic acid = 100 / 1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.03 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.2	1.77	1 - 2	6.00
2	6.7	2.59		