


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

N-FMOC-DL-Phenylalanine

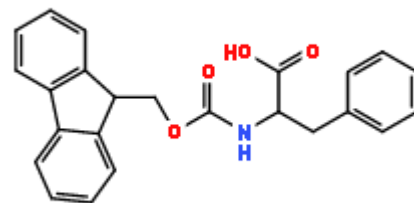

## CHROMATOGRAM

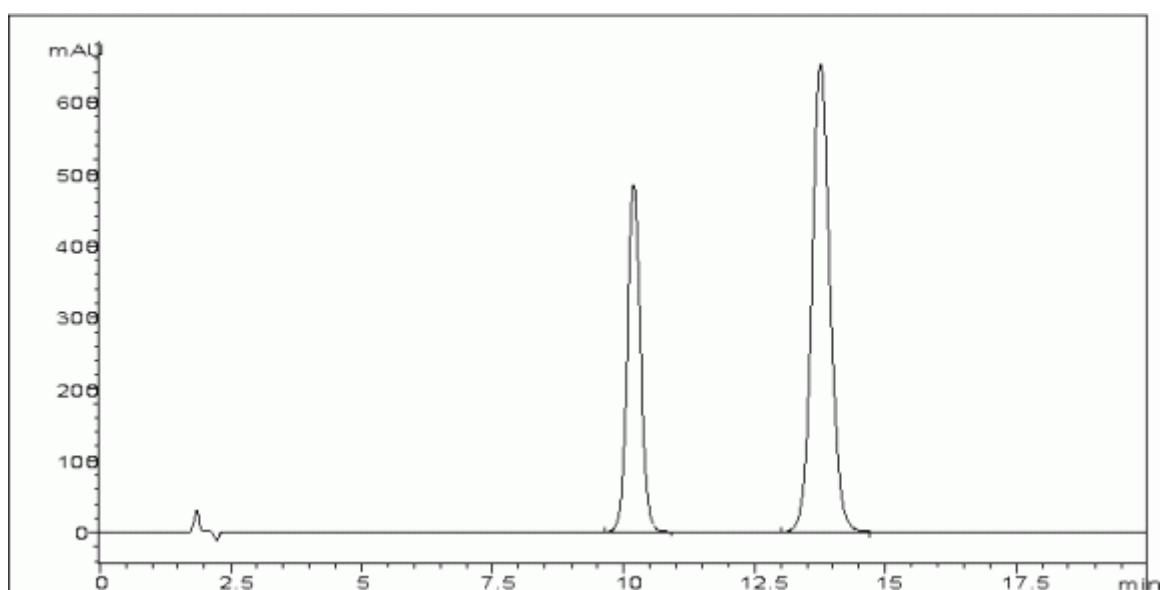

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® QN-AX		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	methanol / acetonitrile / acetic acid = 48.5 / 48.5 / 3		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.2	4.48	1 - 2	6.68
2	13.8	6.40		