


APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name

(±)-Hydrobenzoin

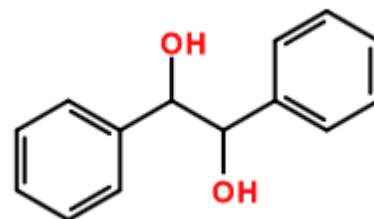

CHROMATOGRAM

NO CHROMATOGRAM

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	CO ₂ / methanol = 80 / 20		
Chromatographic Mode	SFC		
Flow rate (mL/min)	5	Temperature (°C)	35
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	1.8	1.98	1 - 2	1.29
2	2.1	2.57		