


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

1-(Benzyloxycarbonyl)-2-tert-butyl-3-methyl-4-imidazolidinone

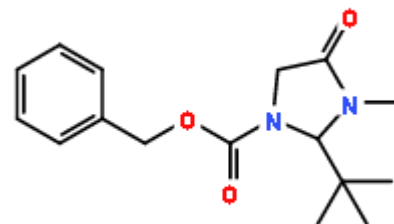

CHROMATOGRAM

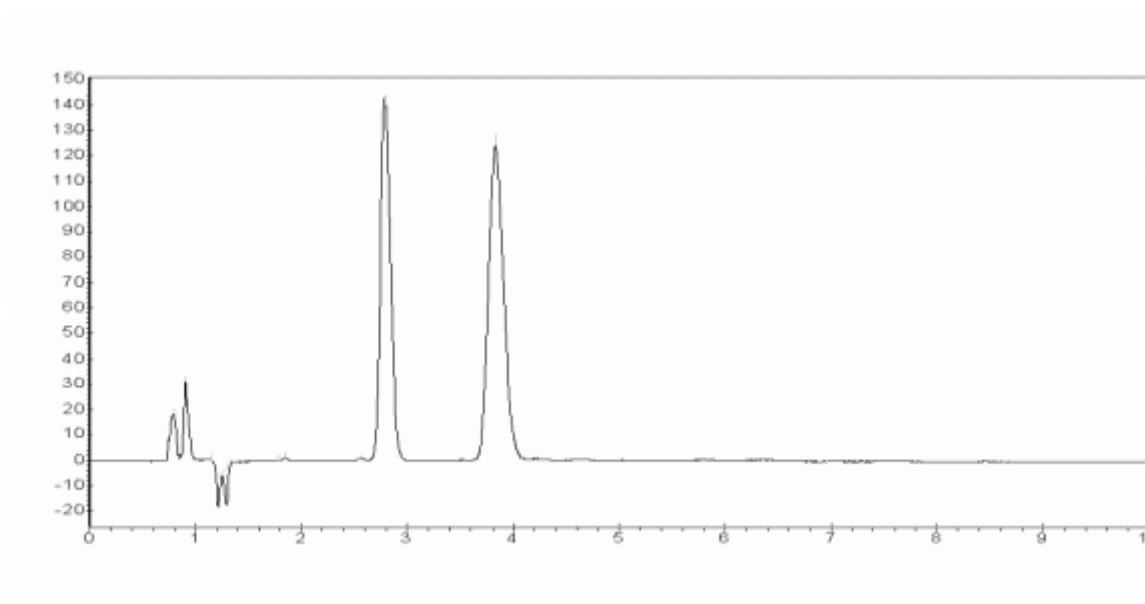

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	CO ₂ / tetrahydrofuran = 75 / 25		
Chromatographic Mode	SFC		
Flow rate (mL/min)	4	Temperature (°C)	35
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.8	2.58	1 - 2	1.52
2	3.8	3.91		