


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Diperodon

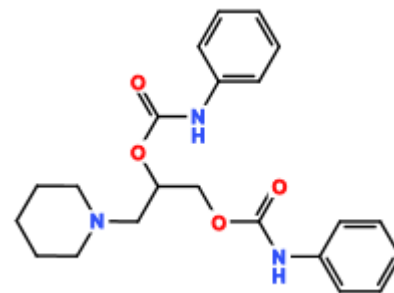

## CHROMATOGRAM

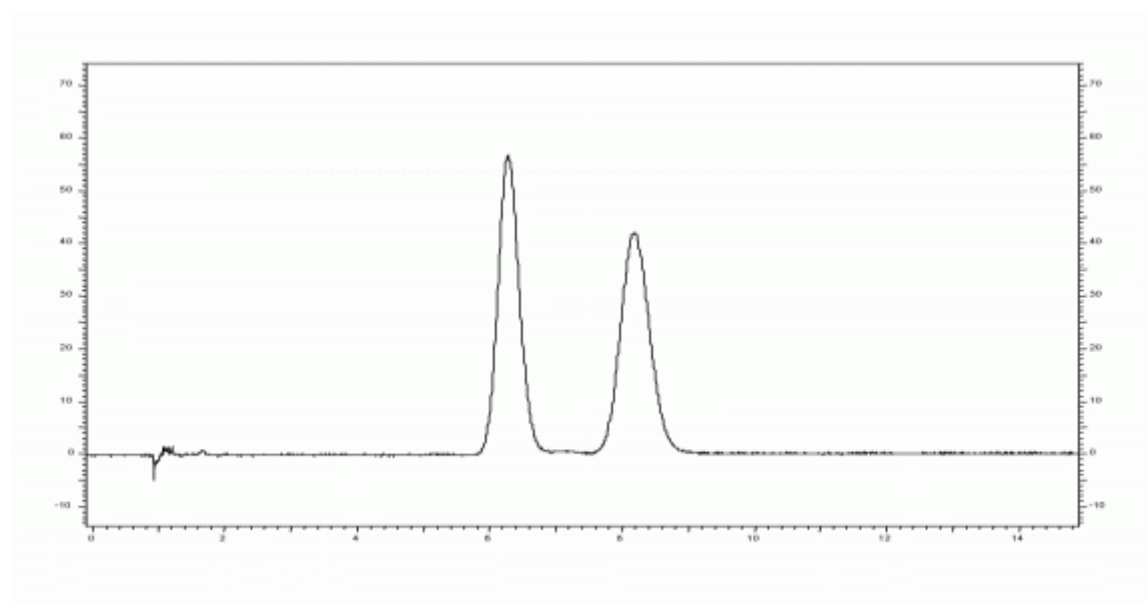

## ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OZ-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	CO <sub>2</sub> / methanol + diethylamine 0.2% = 70 / 30		
Chromatographic Mode	SFC		

Flow rate (mL/min)	3	Temperature (°C)	RT
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.4	5.39	1 - 2	2.65
2	8.3	7.29		