


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Epinephrine

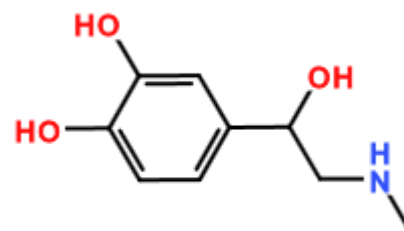

CHROMATOGRAM

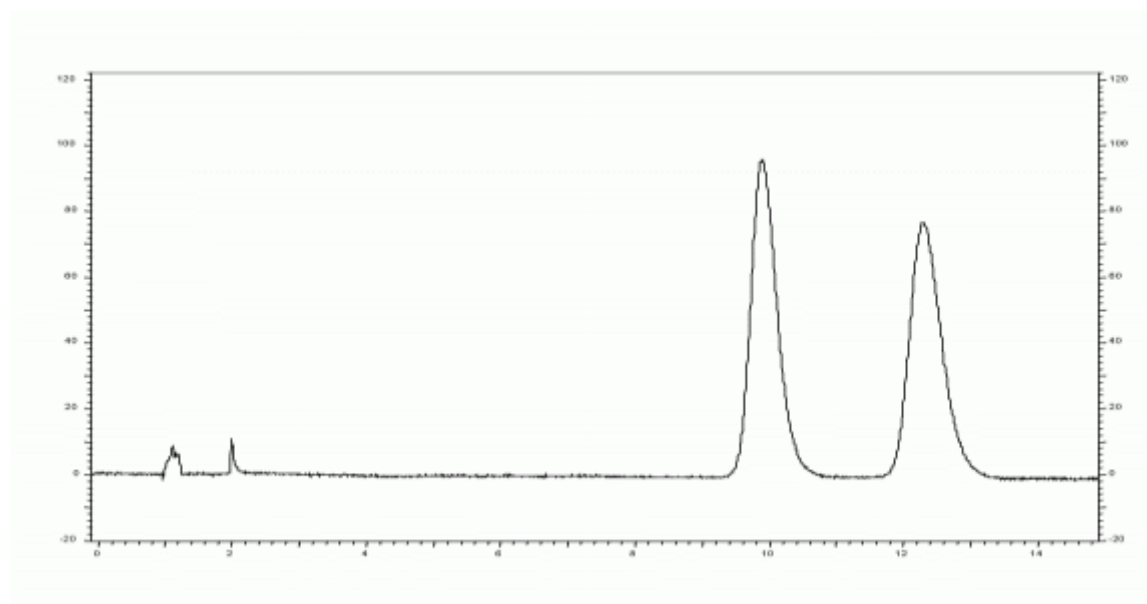

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OZ-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	CO2 / methanol + Ethanesulfonic acid 0.2% = 85 / 15		
Chromatographic Mode	SFC		
Flow rate (mL/min)	3	Temperature (°C)	RT
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.0	9.00	1 - 2	2.97
2	12.4	11.39		