


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

gamma-(4-Fluorophenyl)-gamma-butyrolactone

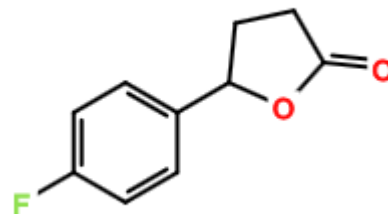

CHROMATOGRAM

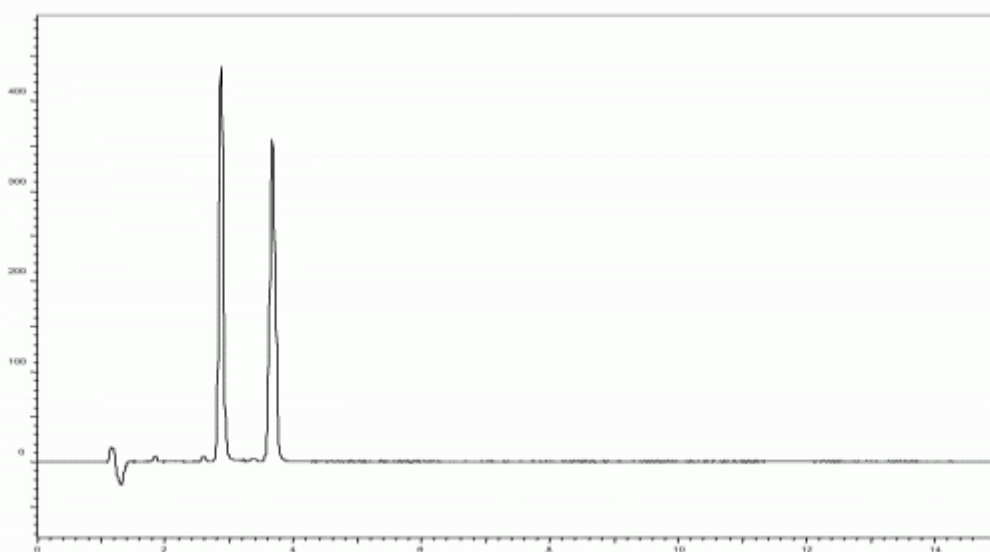

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AY-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	CO ₂ / 2-propanol = 85 / 15		
Chromatographic Mode	SFC		
Flow rate (mL/min)	3	Temperature (°C)	RT
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.9		1 - 2	1.43
2	3.7			5.53