


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Methyl 1-benzyl-5-oxo-3-pyrrolidinecarboxylate

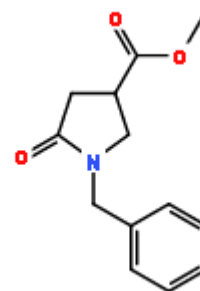

CHROMATOGRAM

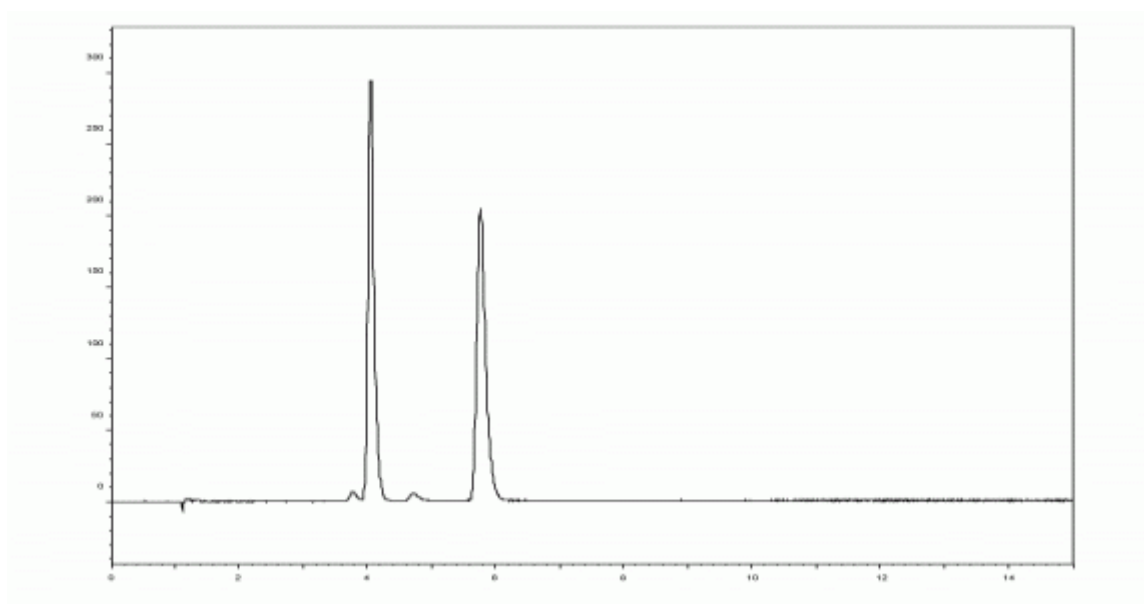

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OZ-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	CO ₂ / methanol = 85 / 15		
Chromatographic Mode	SFC		
Flow rate (mL/min)	3	Temperature (°C)	RT
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.1	3.05	1 - 2	8.20
2	5.8	4.77		