


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

N-FMOC-DL-valine

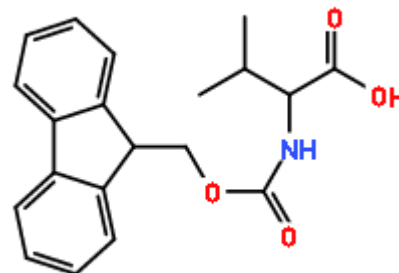

CHROMATOGRAM

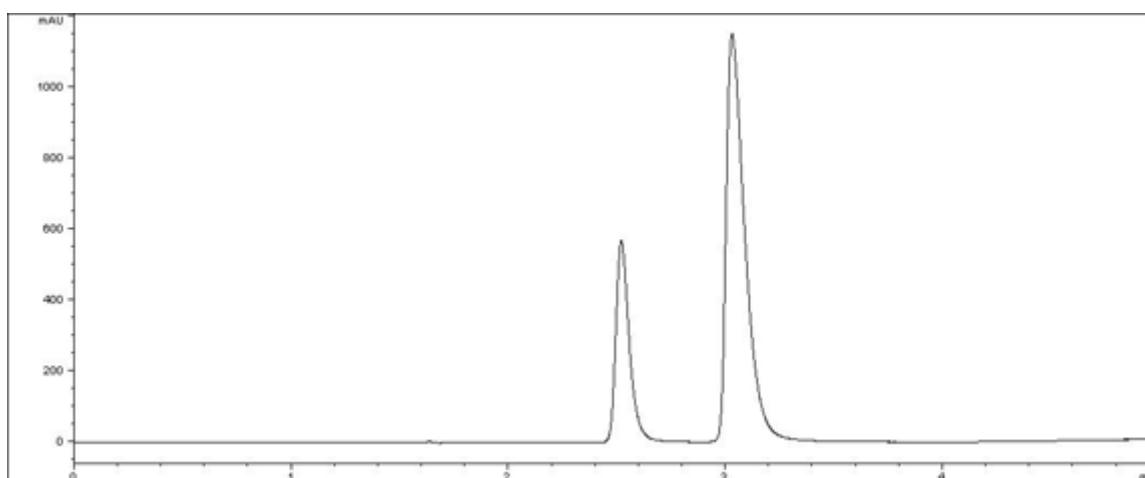

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	25mM acetic acid in methanol / water = 98 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order	D/L		

No.	Rt	K'	Alpha Value	Resolution
1	2.5	0.53	1 - 2	3.98
2	3.0	0.84		