


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Pregabalin

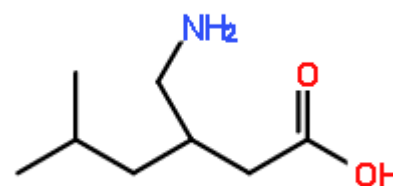

## CHROMATOGRAM

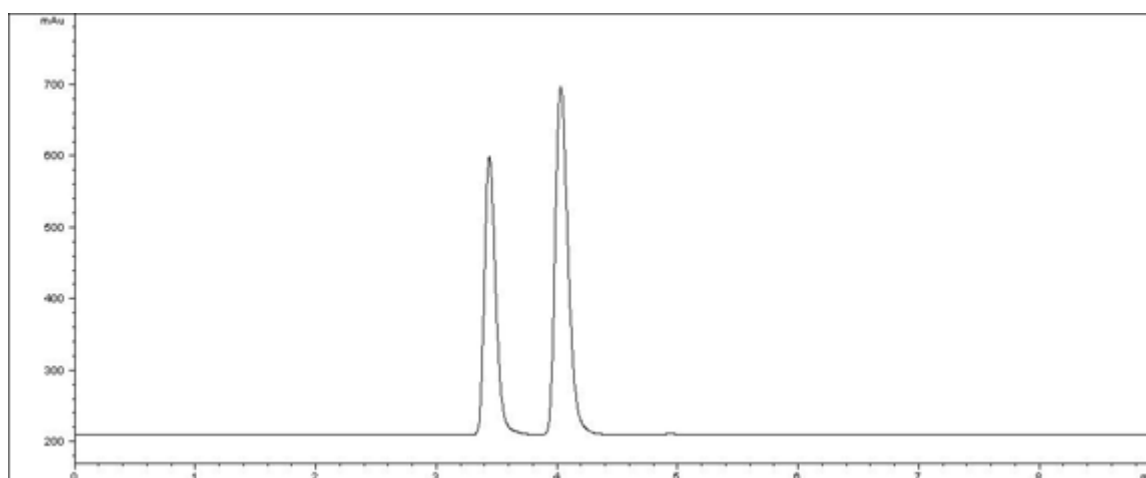

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 90 / 10		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order	R/S		

No.	Rt	K'	Alpha Value	Resolution
1	3.4	0.90	1 - 2	3.32
2	4.0	1.22		