

## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Cyclopenthiazide



## CHROMATOGRAM

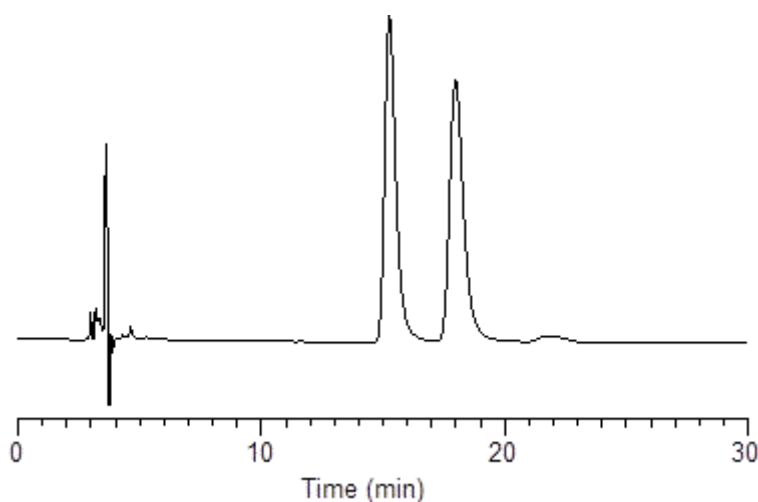

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 80 / 20 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	15.3	4.10	1 - 2	2.99
2	18.0	5.00		