


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Primaquine

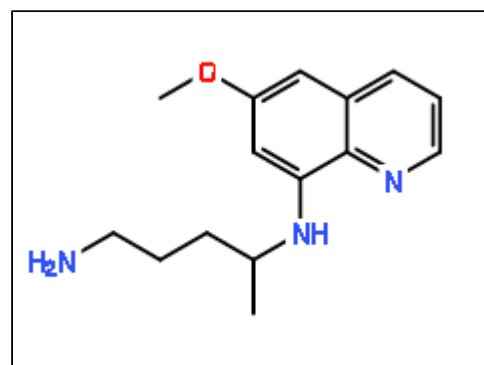

CHROMATOGRAM

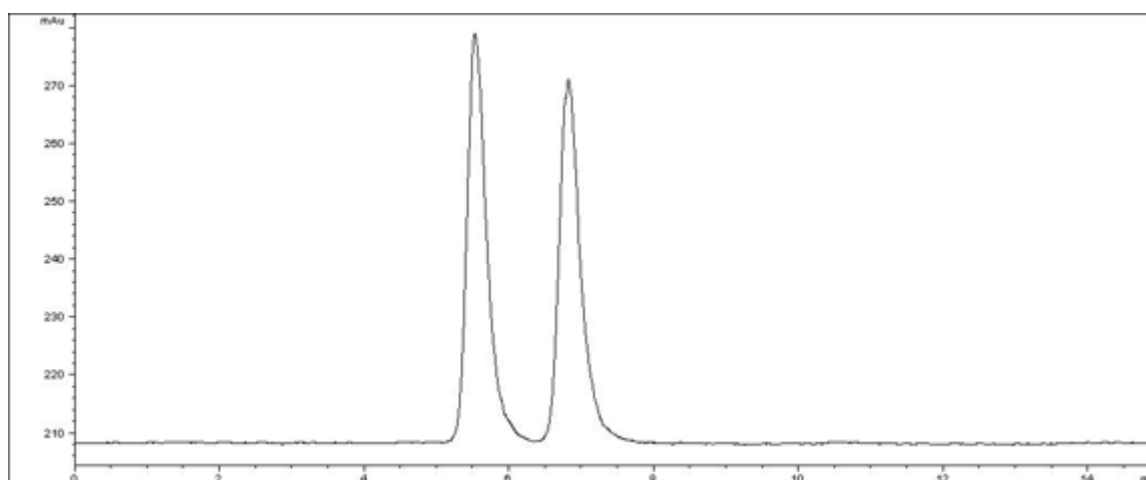

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	12.5mM triethylamine in methanol / water = 98 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0015 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.5	2.00	1 - 2	2.63
2	6.8	2.71		