


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Sulindac methyl ester

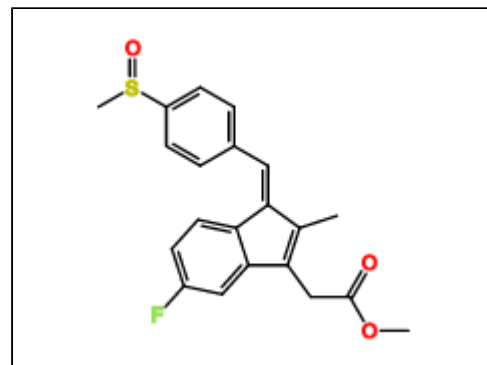

CHROMATOGRAM

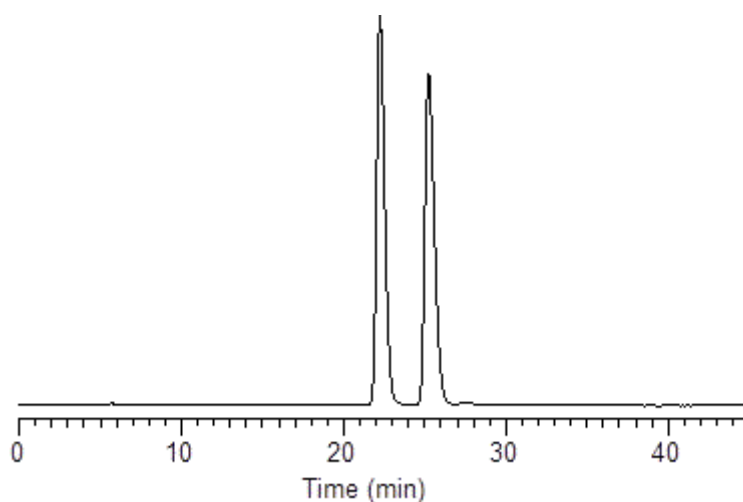

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ID 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 70 / 30		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	22.2	6.40	1 - 2	3.03
2	25.2	7.40		