


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Thioridazine

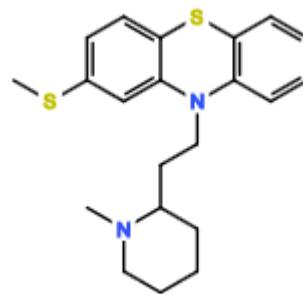

CHROMATOGRAM

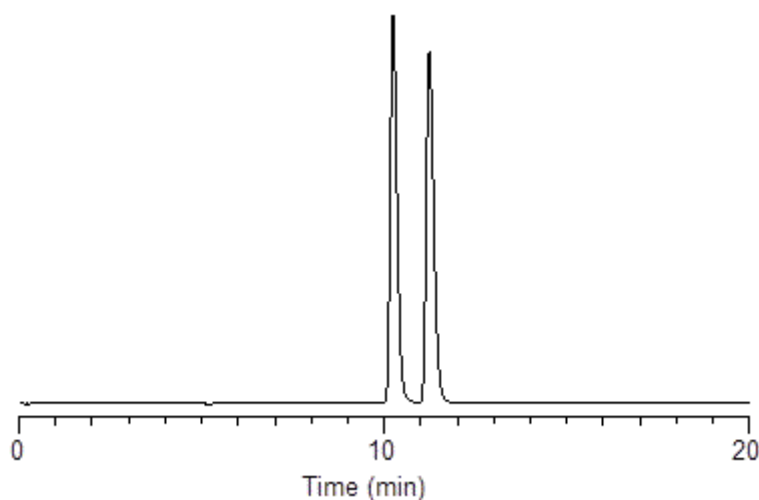

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ID 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / methyl tert-butyl ether / ethanol / diethylamine = 95 / 5 / 1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.003 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.2	2.40	1 - 2 1.14	2.83
2	11.2	2.70		