


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

(±)-N,N'-Bis(alpha-methylbenzyl) sulfamide

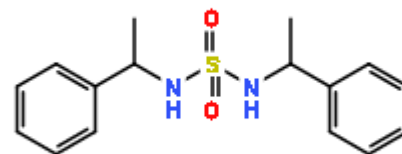

## CHROMATOGRAM

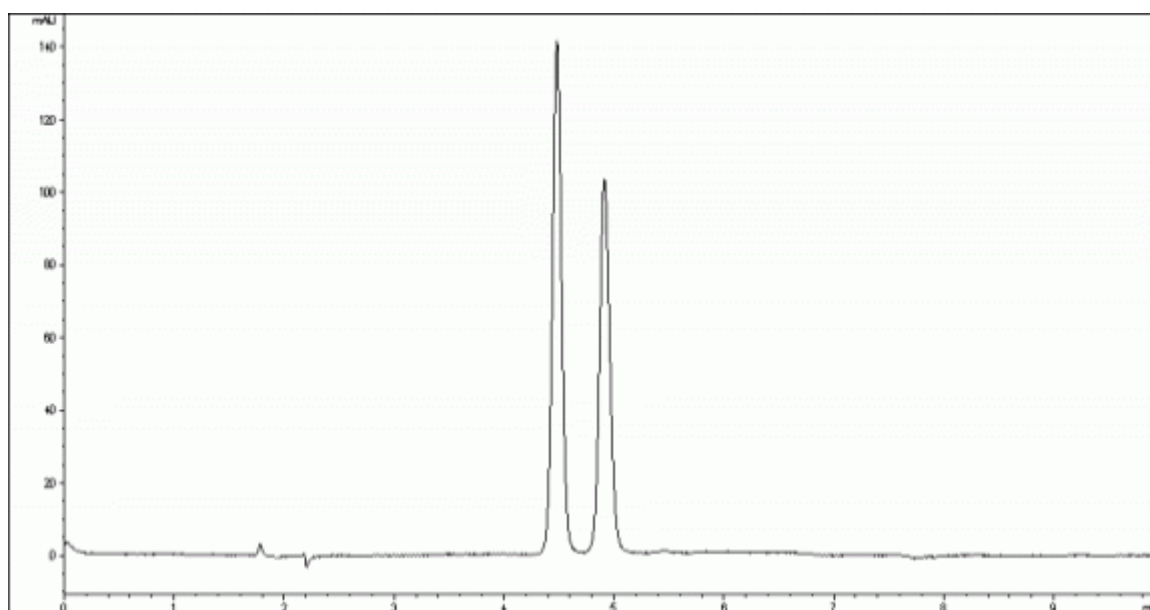

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 80 / 20		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.5	1.63	1 - 2	2.68
2	4.9	1.88		