


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Propranolol

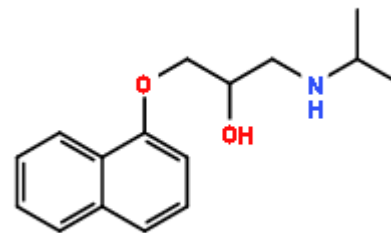

CHROMATOGRAM

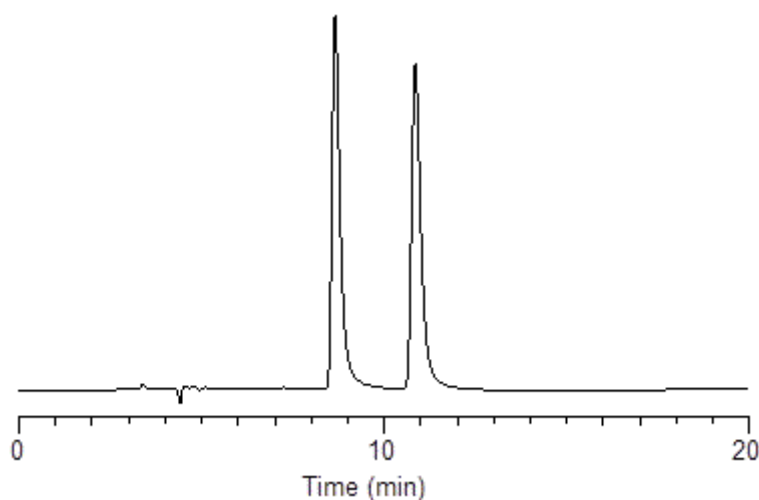

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 93 / 7 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.002 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.7	1.90	1 - 2	5.33
2	10.9	2.60		