


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

(4R,5S/4S,5R)-1,5-Dimethyl-4-phenyl-2-imidazolidinone

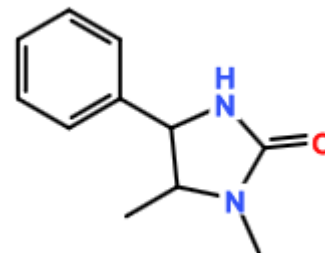

## CHROMATOGRAM

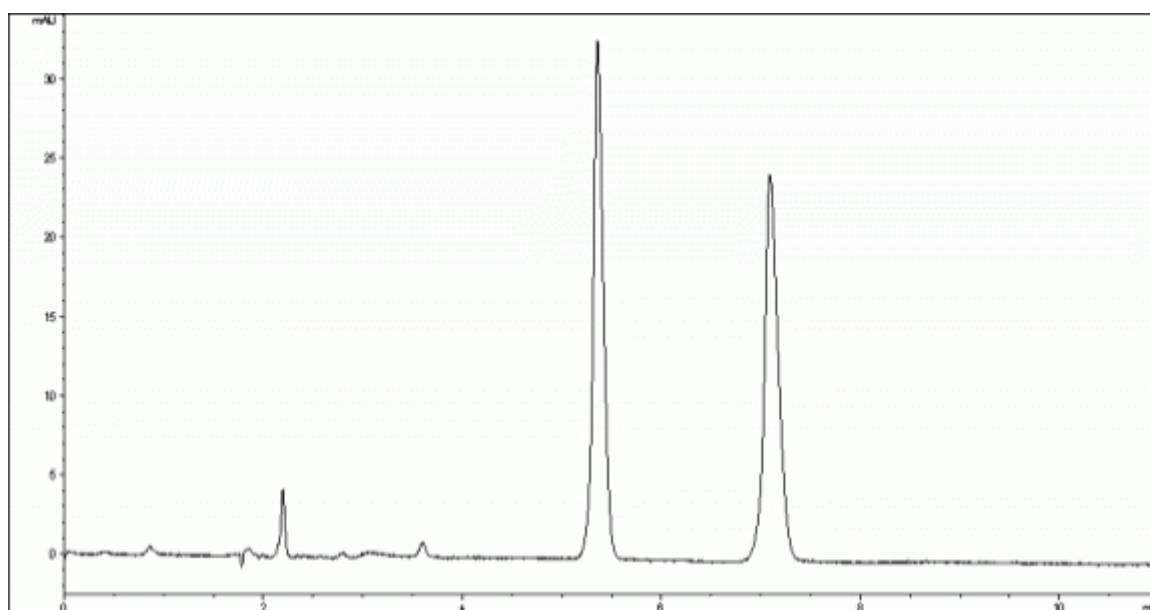

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	dichloromethane / ethanol = 98 / 2		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.4	2.15	1 - 2	8.32
2	7.1	3.17		