


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

1,5-Dihydroxy-1,2,3,4-tetrahydronaphtalene

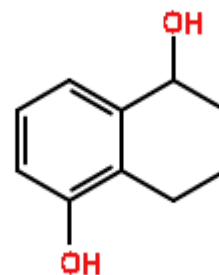

CHROMATOGRAM

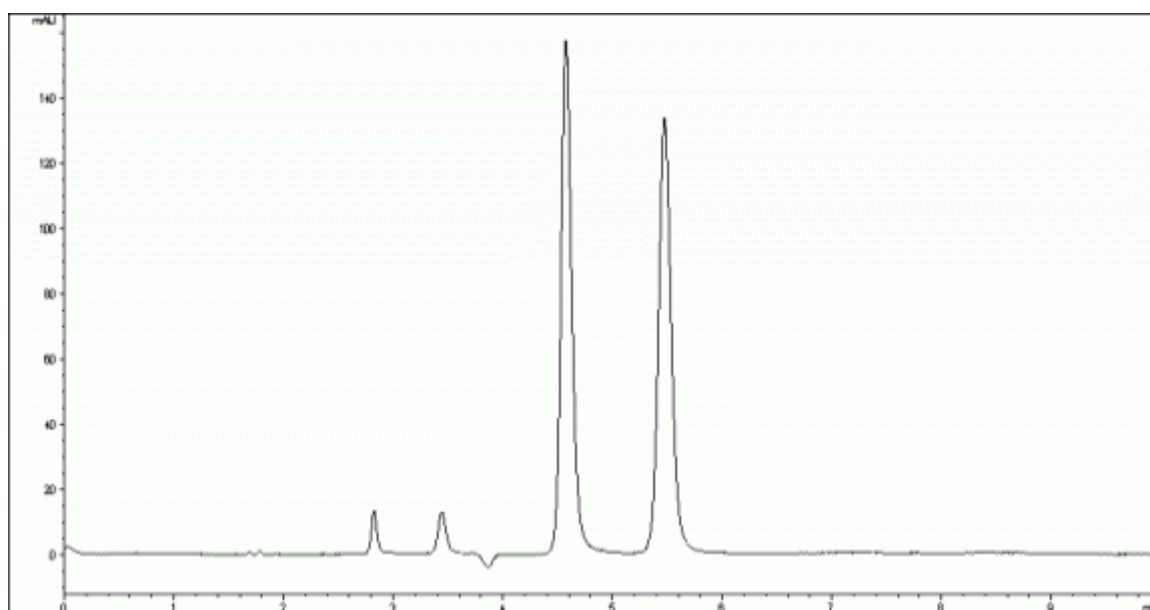

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / ethanalamine = 85 / 15 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.6	1.69	1 - 2	4.78
2	5.5	2.22		