


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Alogliptin

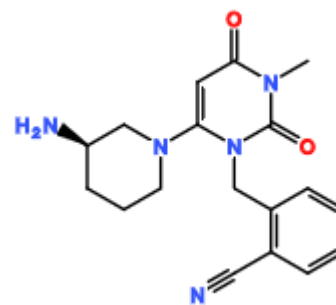

CHROMATOGRAM

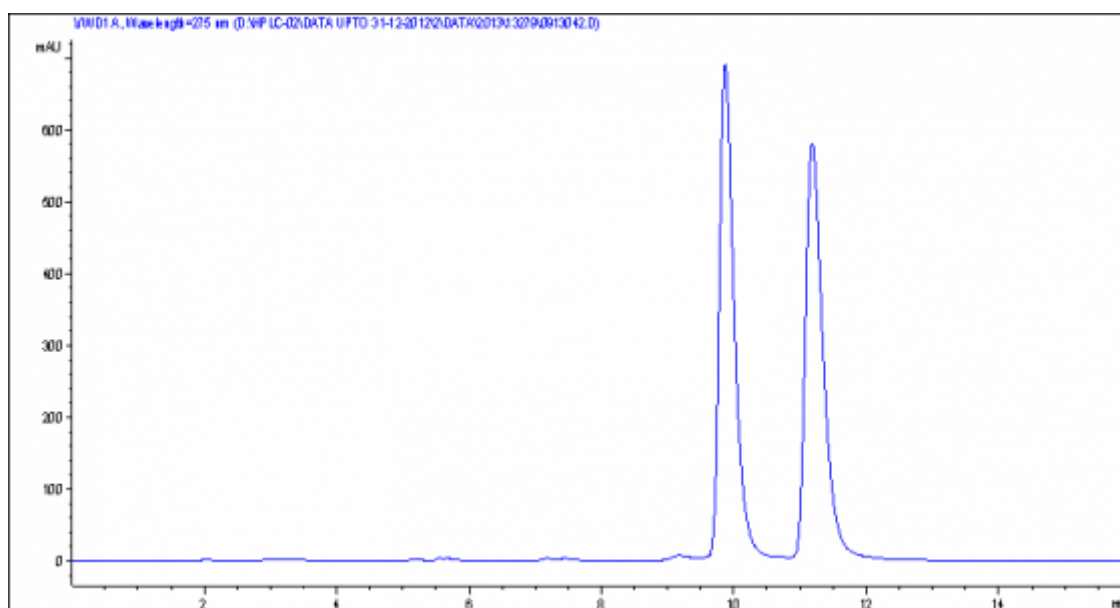

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IF-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	ethanol / methanol / diethylamine = 70 / 30 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.7	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 275 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.9	1.30	1 - 2	3.02
2	11.2	1.61		