


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Apremilast

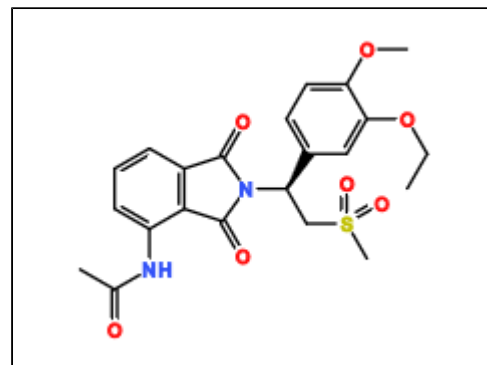

## CHROMATOGRAM

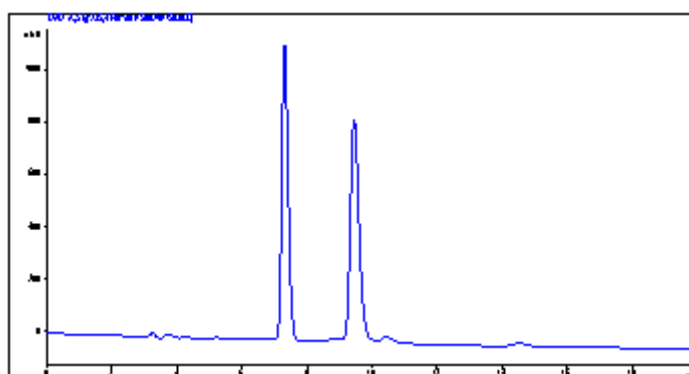

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ID 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methyl tert-butyl ether / ethanol / diethylamine = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 235 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.3	1.44	1 - 2	5.17
2	9.5	2.15		