


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Brinzolamide

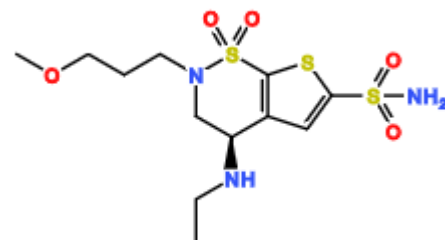

CHROMATOGRAM

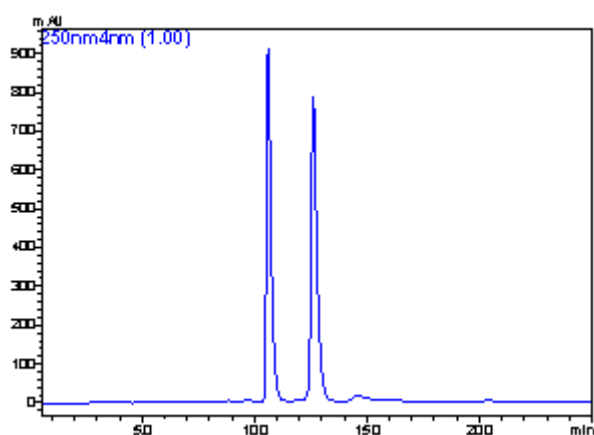

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IE-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 80 / 20 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 250 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.6	2.53	1 - 2	4.94
2	12.6	3.20		