


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Celecoxib

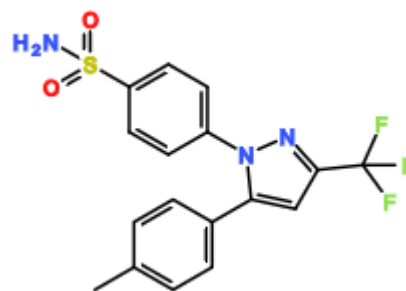

CHROMATOGRAM

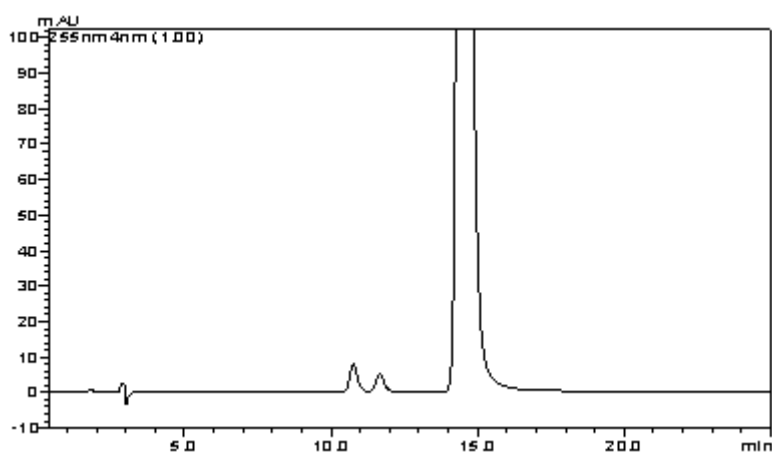

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 92 / 8 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 255 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.8	4.98	1 - 2	1.78
2	11.7	5.47	2 - 3	4.88
3	14.5	7.05		