


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Clodinafop propargyl

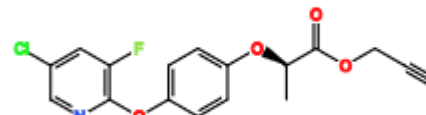

## CHROMATOGRAM

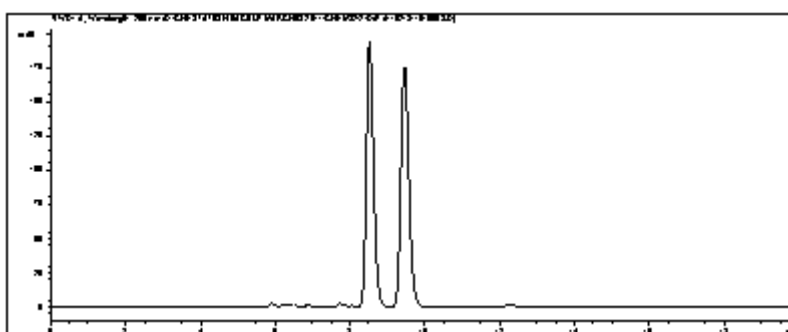

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALCEL <sup>®</sup> OD-H 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / isopropylamine / trifluoroacetic acid = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 265 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.5	1.84	1 - 2	2.54
2	9.5	2.15		