


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Cobicistat

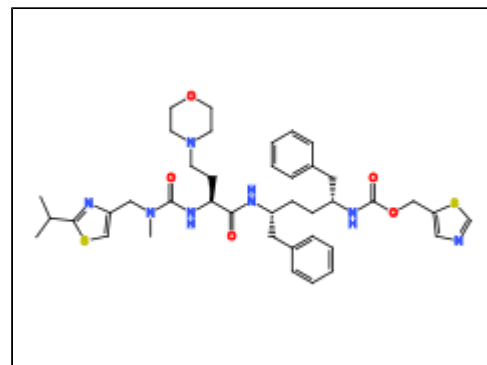

## CHROMATOGRAM

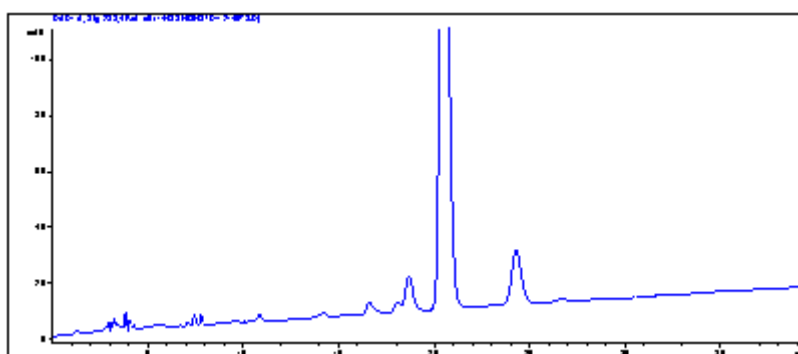

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IE-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 65 / 35 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	40
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 238 nm		
Elution Order	RRR, SSS, API, SRS		

No.	Rt	K'	Alpha Value	Resolution
1	16.6	4.54	1 - 2	2.11
2	18.7	5.23	2 - 3	2.44
3	20.5	5.83	3 - 4	4.31
4	24.3	7.10		