


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Darunavir

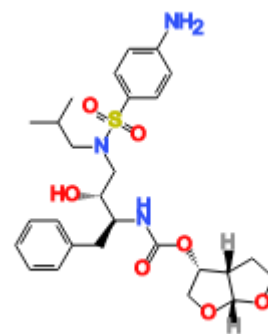

CHROMATOGRAM

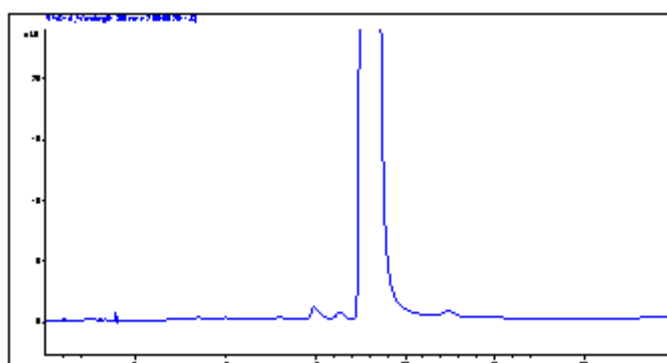

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IE-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 70 / 30 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.03 (mg)	Sample Concentration	
Detection	UV-VIS 265 nm		
Elution Order	RS, SS, API, RR		

No.	Rt	K'	Alpha Value	Resolution
1	14.9	3.98	1 - 2	1.12
2	16.4	4.46	2 - 3	1.11
3	17.8	4.95	3 - 4	1.31
4	22.4	6.47		