


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Eletriptan

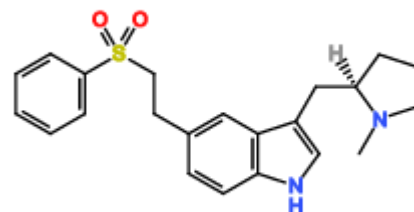

CHROMATOGRAM

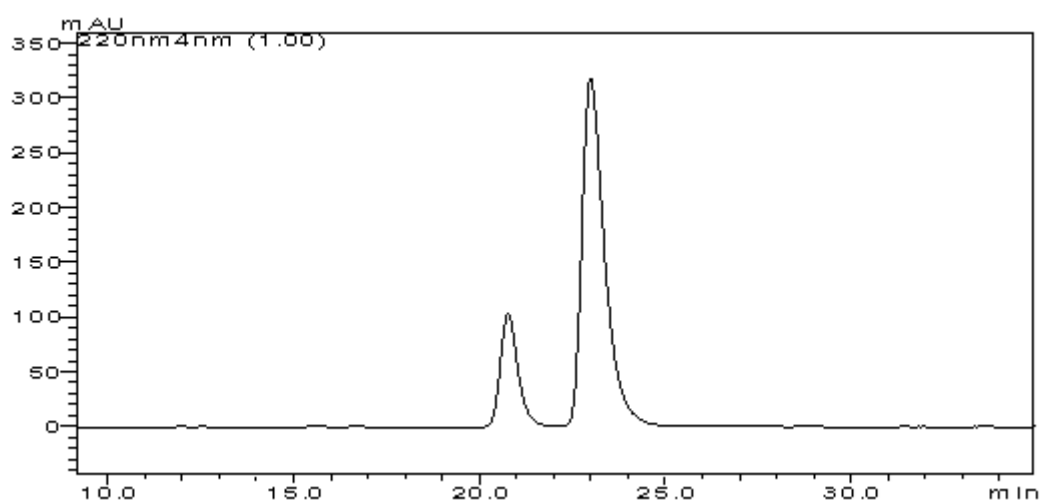

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / methanol / trifluoroacetic acid / triethylamine = 88 / 6 / 6 / 0.3 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	20.8	10.54	1 - 2	1.12
2	23.0	11.77		2.31