


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Linezolid

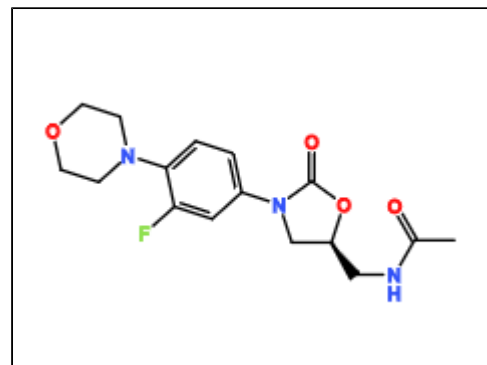

CHROMATOGRAM

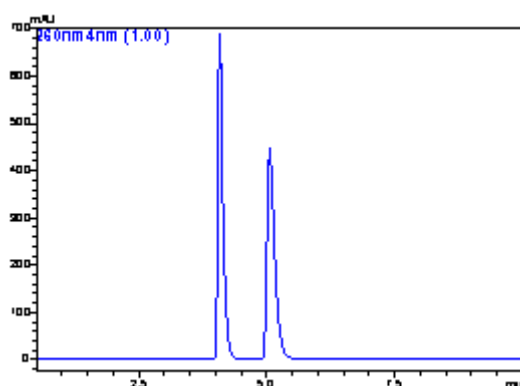

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	acetonitrile / methanol / diethylamine = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	40
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 260 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.1	0.35	1 - 2	3.77
2	5.0	0.68		