


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Nicardipine

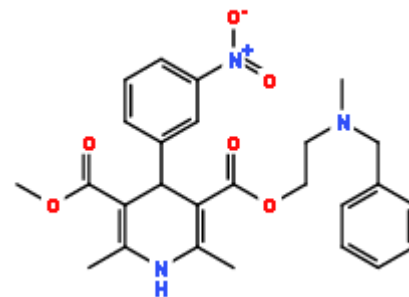

## CHROMATOGRAM

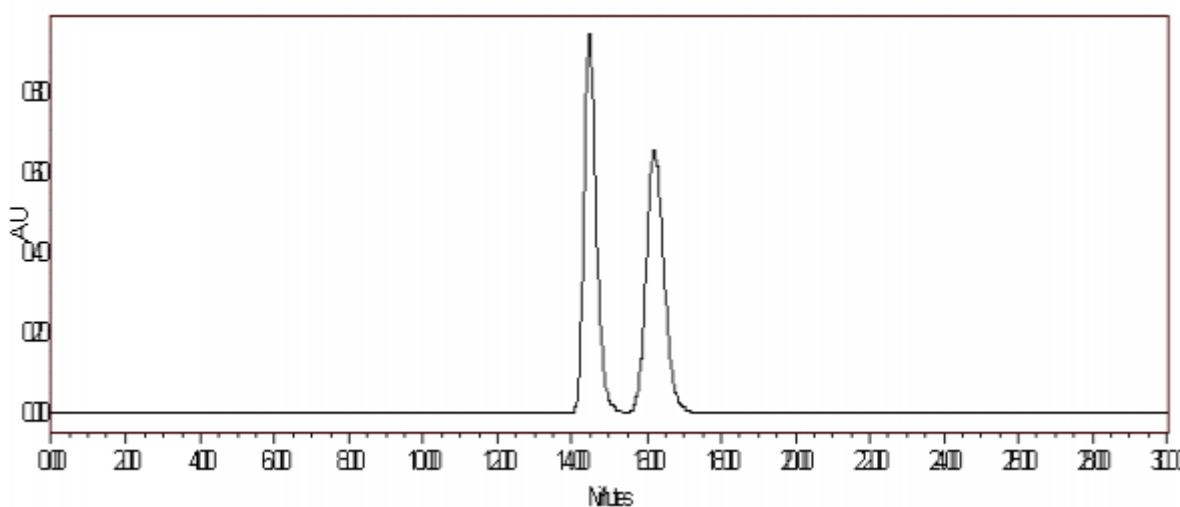

## ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OJ-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	ethanol / acetonitrile / Ammonia = 98 / 2 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	0.3	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 240 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	14.5	0.45	1 - 2	2.30
2	16.2	0.62		