


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Penicillamine

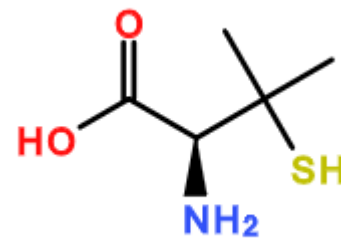

CHROMATOGRAM

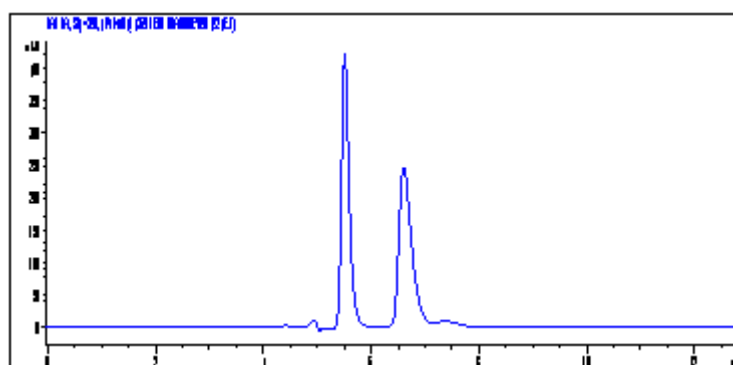

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CROWNPAK [®] CR-I(+) 3 x 150 mm (5 μm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 80 / 20		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.2	Temperature (°C)	10
Injection Volume (μL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.5	0.28	1 - 2	3.27
2	6.6	0.54		