


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Pralatrexate

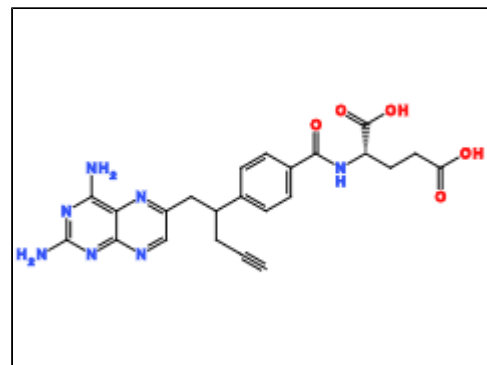

CHROMATOGRAM

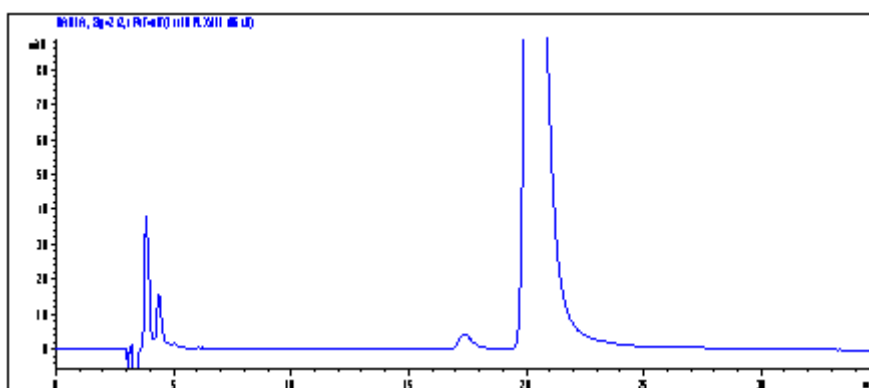

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ID		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / methanol / trifluoroacetic acid / n-butylamine = 80 / 10 / 10 / 0.5 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	40
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 242 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	17.5	4.83	1 - 2	1.20
2	20.4	5.80		2.82