


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Rasagiline

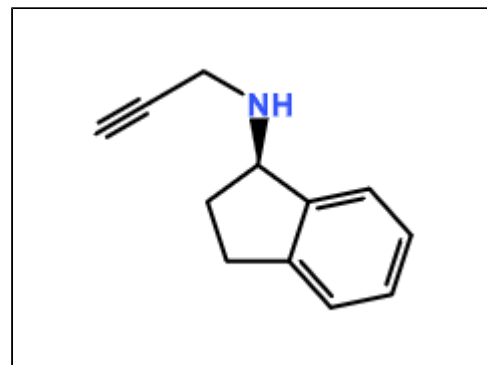

## CHROMATOGRAM

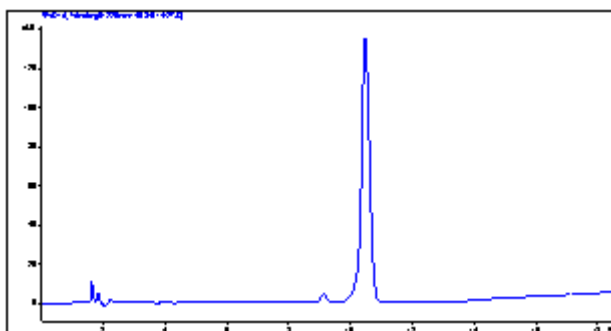

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IF-3 4.6 x 150 mm (3 µm)		
Mobile Phase Composition	10mM ammonium bicarbonate aq. / methanol / acetonitrile = 50 / 40 / 10		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.1	4.08	1 - 2	3.21
2	10.5	4.82		