


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Silodosin

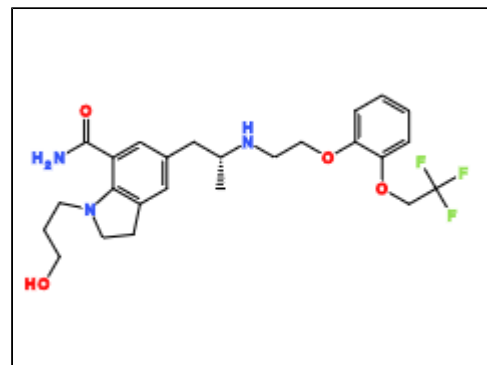

CHROMATOGRAM

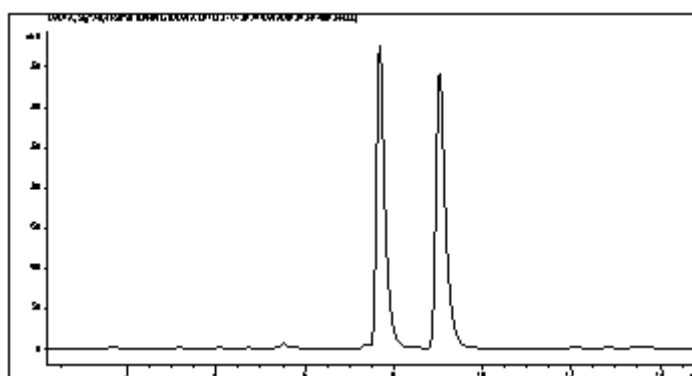

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IF-3 4.6 x 250 mm (3 µm)		
Mobile Phase Composition	methyl tert-butyl ether / ethanol / diethylamine = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.7	1.56	1 - 2	3.87
2	9.0	2.01		