


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Valacyclovir

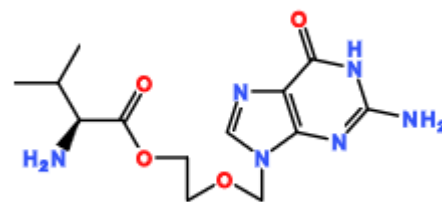

CHROMATOGRAM

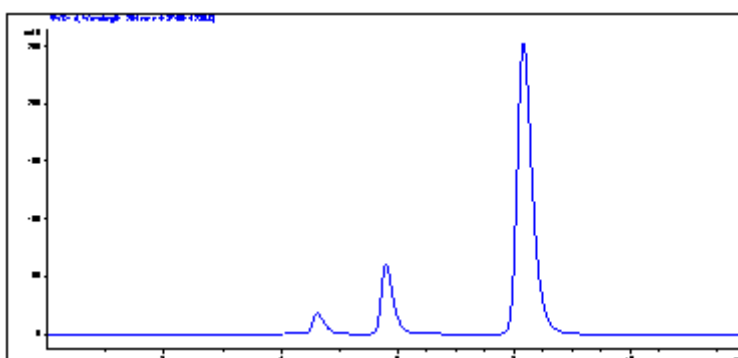

ANALYTICAL REPORT

Column Trade Name	CROWNPAK® CR-I(+)		
Column Size (Particle Size)	3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 80 / 20		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.2	Temperature (°C)	25
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order	D, L		

No.	Rt	K'	Alpha Value	Resolution
1	5.8	0.35	1 - 2	5.64
2	8.2	0.90		