


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

2,3-Dibenzoyl-DL-tartaric acid

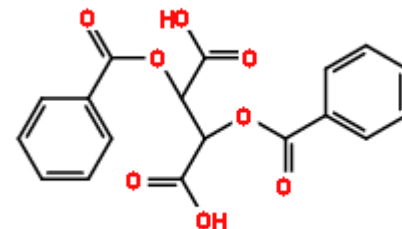

CHROMATOGRAM

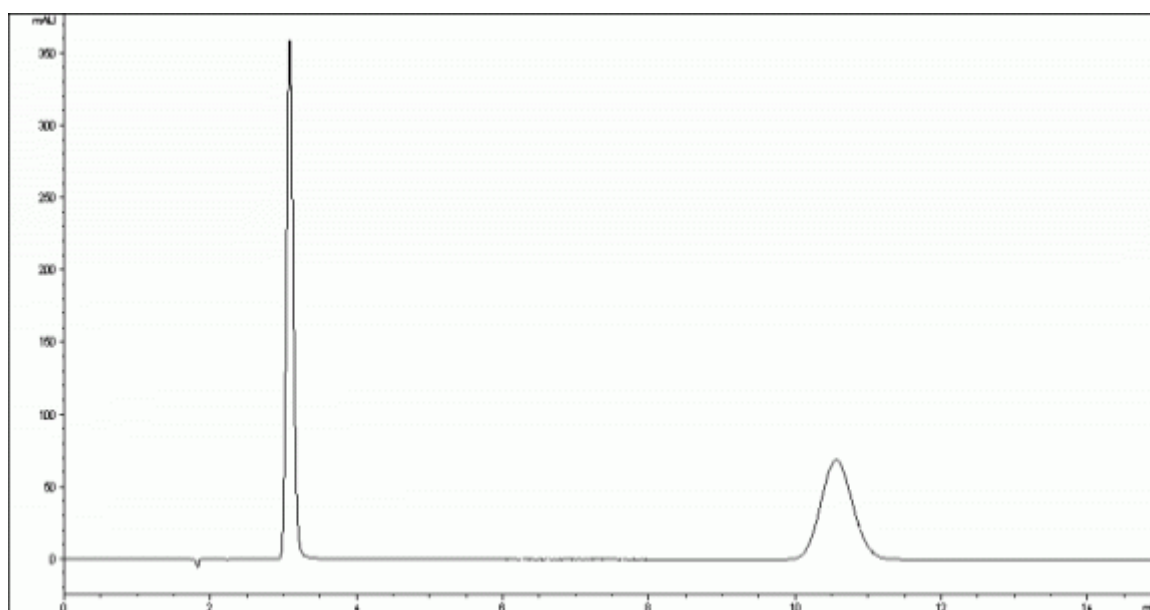

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	dichloromethane / ethanol / trifluoroacetic acid = 98 / 2 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.1	0.80	1 - 2	6.44
2	10.6	5.18		
				15.27