


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Cinacalcet

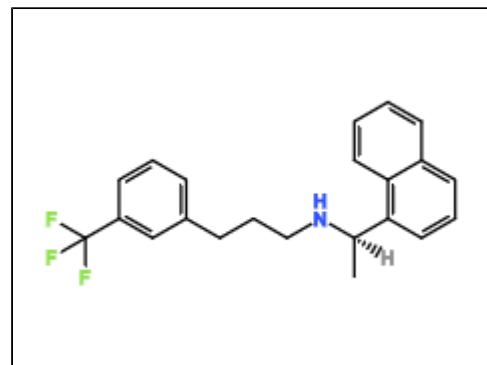

CHROMATOGRAM

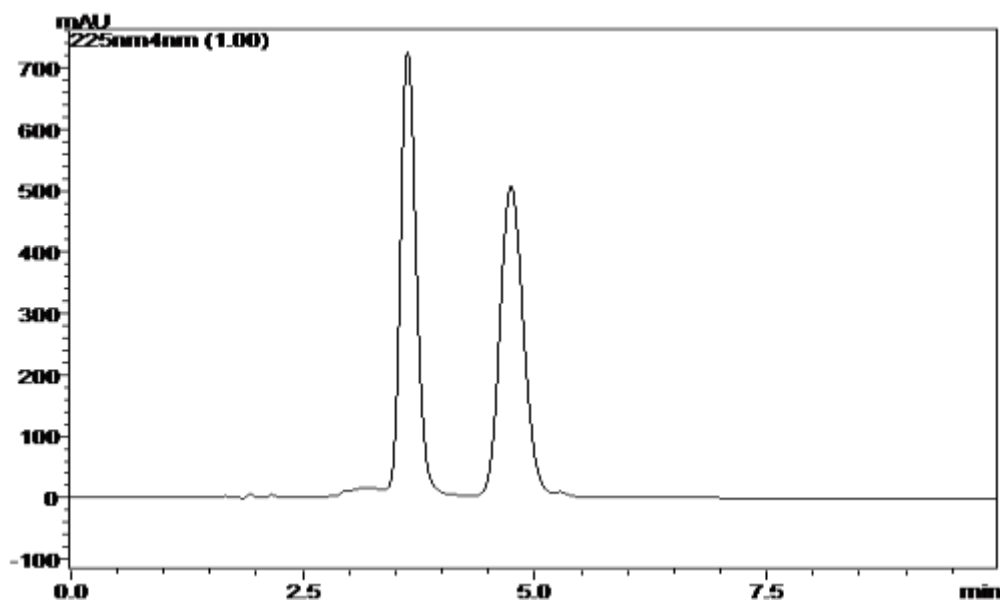

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-3R		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	50mM ammonium bicarbonate aq. / methanol = 10 / 90		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 225 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.6	1.02	1 - 2	2.86
2	4.7	1.63		