


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Fluvastatin sodium

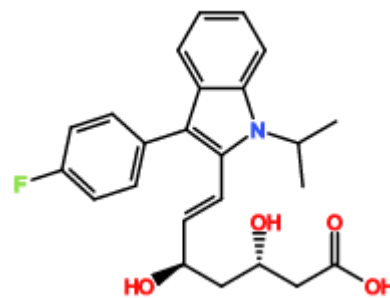

CHROMATOGRAM

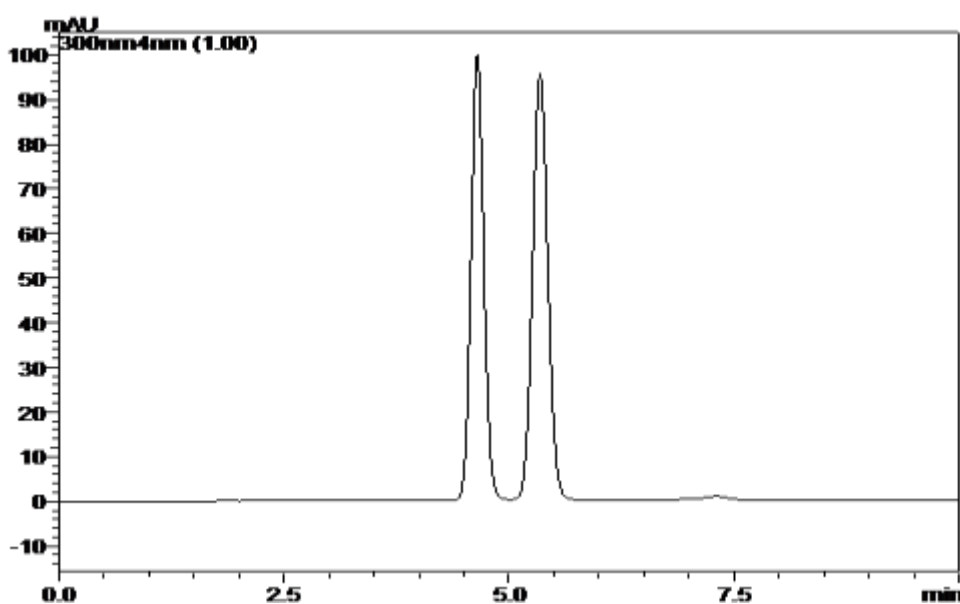

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OZ-RH		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	0.1% Formic acid aq. / acetonitrile = 55 / 45		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 300 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.6	1.58	1 - 2	2.34
2	5.3	1.97		