


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Garenoxacin mesylate

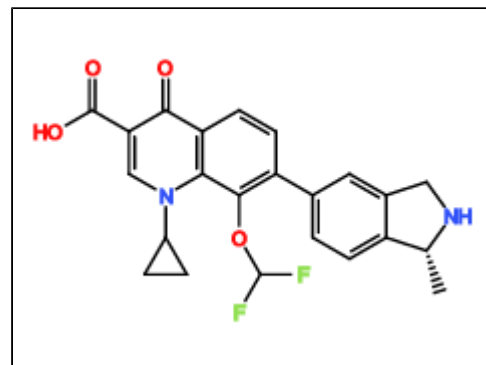

CHROMATOGRAM

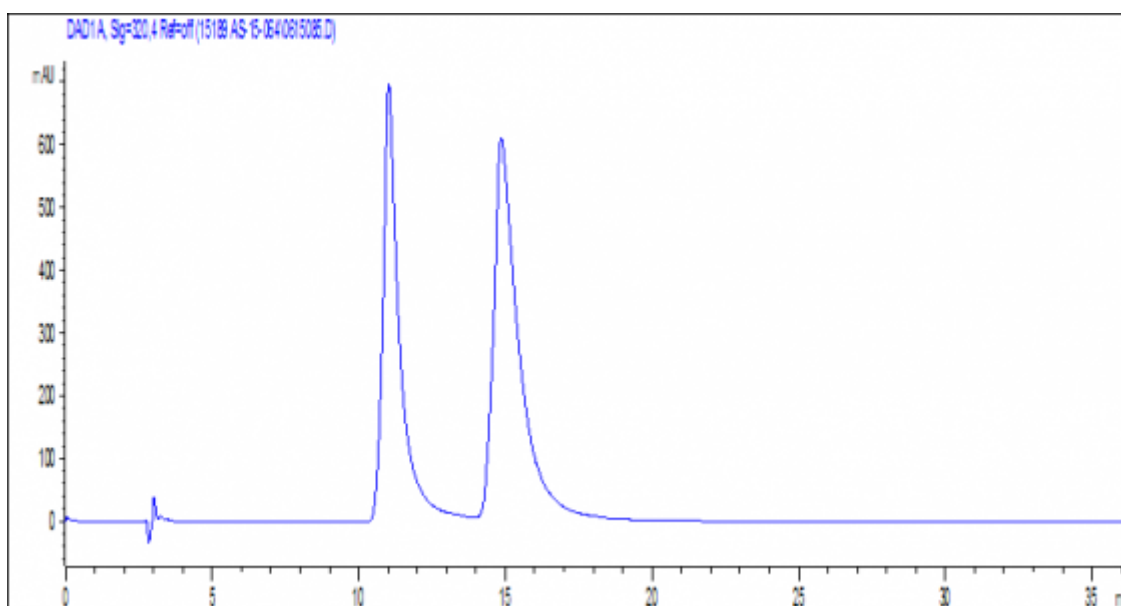

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	100mM formic acid + 25mM diethylamine in methanol / tetrahydrofuran / water / acetic acid = 49 / 49 / 2 / 0.5		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	15
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 320 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.0	2.06	1 - 2	1.52
2	14.9	3.13		3.50