


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

4-Borono-DL-phenylalanine

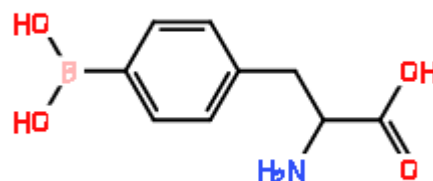

## CHROMATOGRAM

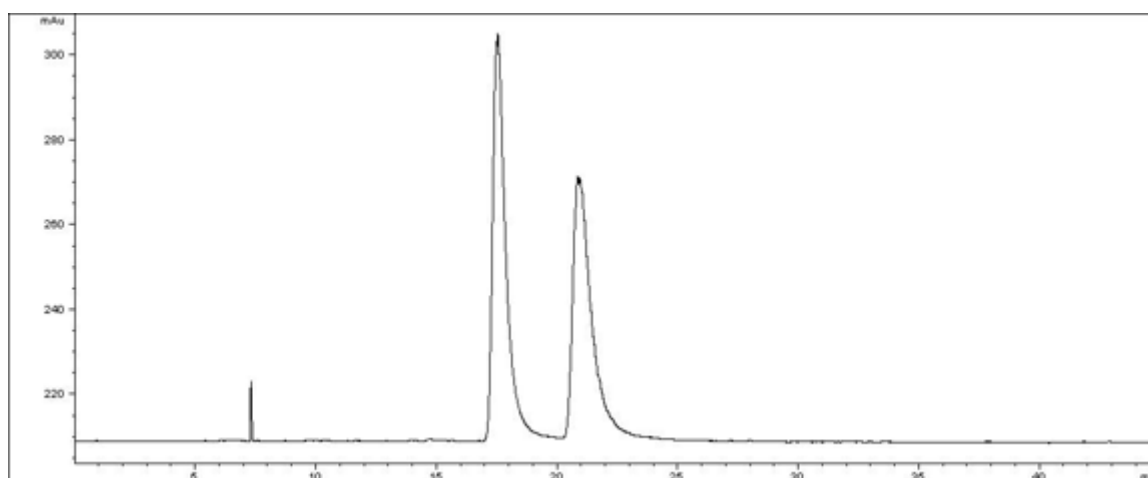

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	4 x 250 mm (3 μm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / acetonitrile / water = 49 / 49 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	10
Injection Volume (μL)	0.005 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	17.5	2.62	1 - 2	2.98
2	20.8	3.30		